

Coaxmed

General realties	
Main voltage	240 Vac 50/60Hz
Maximum Power	450 VA
Fuses	2 x F4A 250 V
Maximum Vacuum	66.6kPa
Minimum Temperature Cryo	-10°
Safety Class and Type Applied Part	I BF
Risk Classification Based on 93/42/CEE	Class IIb
Dimensions [W x D x H]	680x470x1400mm
Maximum weight with handpieces	63kg

Frequency	28kHz
Maximum Power	35W
Operating Mode	Continuous and pulsed
Activation	Footswitch


N° Diode Laser	4
Wavelenght	808nm ± 5
CW output power CW for Diode	1W

Maximum Output Voltage	120Vpp @ 1kohm
Maximum Current Density	1.5mA/cm²
Single Pulse Duration	< 1ms

Function	Bipolar/Monopolar
Frequency	From 300 kHz to 1000kHz
Waveform	Amplitude Modulated Signal
Maximum Power	100W
Operating Mode	Continuous and pulsed
Activation	Footswitch

Function	Bipolar/Monopolar
Frequency	From 300 kHz to 1000kHz
Waveform	Amplitude Modulated Signal
Maximum Power	130W
Operating Mode	Continuous and pulsed
Activation	Footswitch



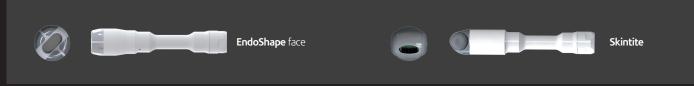


Biotec Italia guarantees the genuineness of a truly Italian product and the superior quality of its products for the infinite attention paid to every single detail of the manufacturing process.

- Biotec Italia Guarantee:
 Entirely made in Italy.
 Made using only Italian semi-finished products.
- Made with top quality materials.
- Made adopting Italian technological expertise and know how.
 Made using Biotec Italia's own technological designs.



COOXMED

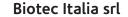

Viale della Repubblica, 20 - 36031 Dueville (Vicenza) tel. 0444.591683 - fax. 0444 361032













Cooltec®: Cools the handpiece's tip preserving the skin surface from overheating. + 3Density®: Automatically modulates field frequency **Flexi®** to vary the depth of penetration of thermal energy in Experience the benefits.Improve results. the skin, from 300KHz to 1MHz. (+) Propatch: single-use, disposable, saline gel patch that The synergistic effect of different technologies that guarantees uniform thermal emission, eliminating the effectively treat face and body aesthetic needs, as risk of superficial blisters. non-surgical body contouring, cellulite reduction, skin tightening, and wrinkle reduction. + EndoShape®: cells stimulation by adjustable rhythmic suction that facilitates R F energy penetration Combination of rhythmic vacuum, microneedling and radiofrequency to induce dermal coagulation Vaculase[®] **Flexi**® Phase1: Deep Laser Stimulation Heats up the deeper dermis at 360 degrees around each microneedle tip generating collagen tightening The laser diodes reactivate and boosts the arterial, venous and lymphatic micro-circulation, encouraging + Sterile disposable multi-needles electrode from 0.25mm neoangiogenesis of the intradermal to 2.5 mm length. vessels. The laser light energy also has beneficial effects on fibroblast activity. **DuoShape®** Phase2: Localised Cooling 3Density®: Automatically modulates field frequency to The localised cooling system boosts vary the depth of penetration of thermal energy in the the volume-reducing effect, skin, from 200KHz to 600KHz. optimising reabsorption of the edema by stimulating vascular workout. Exclusive simultaneous emission of high-intensity radiofrequency and high-power 28 KHz ultrasound + Phase3: Lipomassage generator which damages fat cell membranes thanks The rhythmic aspiration massage to the cavitational effect. stimulates the lympho-drainage and tissue stimulation processes, PerfectTune® technology cools the ultrasound enhancing restoration of the skin's transducer with active water flow which normal collagen elastic tension. preserves stability and guarantee efficency. Reduction of adipose tissue from 20% to 30% following each session. (+) Immediate return to the normal activities. Possibility to treat different areas of adiposities in the same session. Cryoliposculpt® # High safety level during the treatment. (+) Elimination of the "Ice Block" effect in the treated tissue area. # Restores skin's elasticity and tissue trophism. Reduced discomfort during the treatment.



maximize the investment while minimizing cost

- Coaxmed is a unique technological innovation combining several technologies, the key to an amazing increase in the quality of the aesthetic results.
- Obsolescence is one of the major problems in technology: The need to purchase the latest piece of equipment is the recurring complaint of physicians worldwide.
- Coaxmed is designes to evolve. An upgradable platform with the ability to be upgraded to the latest innovations and technologies in the market without the need of purchasing a new model.

Treatment Indications

- **Stretch** marks and scars
- **4** Alopecia
- **Example 2** Superficial and deep wrinkles
- Acne and open pores

- **Skin laxity**
- ① Pefs (cellulites)
- **Excess fat**
- **Hyperhidrosis**













oShape) Before



e) Before

30 days after 4 sessions

Tissue Active® is the new fat freezing procedures.

3Density[®] is the complex algorithm to deliver more homogeneous heating to multiple tissue depths.

VacuPulse® e Coolfrequency® The innovative way to contour the body by freezing the fat and thightening the skin.